

**THE  FRANKENMAN**  
CS STANDARD PORTFOLIO



 **FRANKENMAN**



WHEN YOU PUT PATIENTS' HEALTH FIRST,  
**REMARKABLE THINGS HAPPEN**





are very proud of our Circular stapler portfolio. Simply the most extensive of any manufacturer on the market.






Many devices have been designed with the input of Global Key opinion leaders to bring unique iterations and innovative devices to the surgical community whilst providing improved functionality, ergonomics and the continual evolution of our circular stapler portfolio, embodies our ethos of "When you put patients' health first, remarkable things happen"



### CS STANDARD

-  420mm Standard length.
-  21, 25, 28 & 32mm.

### SAFETY & EFFICACY




-  Longer, 5.0mm, open staple leg height, to ensure correct "B" Shaped formation & security,
-  even in thicker tissue.
-  Unique engineering of integrated staple driver & cutting blade, gives greater stability & strength leading to improved performance during the firing cycle.
-  Robust Safety locking mechanism prevents premature firing.
-  Clear indicators for opening & closing.

Extended pursestring tying area, reduces risk of inaccurate application.

Reduced "force to fire" vs alternative devices, due to design of washer in anvil, reduces distal movement, & increases stability during the firing cycle. This may also reduce the risk of incomplete firing.






## ANVIL HOLDER

-  All Frankenman Circular Staplers are supplied with an anvil holding rod as standard.
-  Allows secure handling of the anvil when applying a pursestring suture.
-  Controlled handover and docking of the anvil during TransAnal procedures.



## VARIABLE TISSUE COMPRESSION

-  Ensures appropriate staple formation for a variety of tissue thicknesses.
-  "White", "Blue" & "Green" staple closure heights in one device.
-  Gives the surgeon greater control.



2.0MM

1.0MM



# ORDERING CODES AND SPECIFICATIONS

PRODUCT CODE	DESCRIPTION		TOTAL WORKING LENGTH	OUTER DIAMETER (mm)	LUMEN SIZE (mm)	STAPLE QUANTITY	OPEN STAPLE HEIGHT	CLOSED STAPLE HEIGHT	STAPLE CROWN	PER BOX
CS21L	LONG ENDOSTEER™	●	520MM	21	13	16	4.5MM	1.0 - 2.0MM	4MM	3
CS25L	LONG ENDOSTEER™	○		25	17	20	5.0MM			
CS28L	LONG ENDOSTEER™	●		28	20	24				
CS32L	LONG ENDOSTEER™	●		32	24	28				
CS21	STANDARD	●	420MM	21	13	16	4.5MM			
CS25	STANDARD	○		25	17	20	5.0MM			
CS28	STANDARD	●		28	20	24				
CS32	STANDARD	●		32	24	28				

## ALSO AVAILABLE FROM FRANKENMAN

STANDARD CIRCULAR STAPLERS



COMPACT CIRCULAR STAPLERS



LINEAR STAPLERS



LINEAR CUTTING STAPLERS



TRANSCUT™ - TRANSVERSE LINEAR CUTTING STAPLERS




CIRCULAR STAPLERS FOR RECTAL PROLAPSE & ANOPEXY



UNIVERSAL TRANSANAL ACCESS PLATFORM



ENDOSCOPIC LINEAR CUTTER

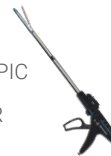

ACCESSORIES



[frankenman.com/academy](http://frankenman.com/academy)

Frankenman are proud to offer several training courses via the Frankenman Academy, in numerous surgical disciplines, with world renowned expert faculty.


These courses are staged to cover both didactic elements and "hands on" skills workshops.

## QUALITY CONTROL

We pride ourselves on world leading continual quality improvement processes.

ALL stapling devices and reloads have a unique serial number, as well a batch/lot number, to allow for more accurate tracking and traceability.

In addition we test fire 100% of our products (before they are then reloaded and packaged), to verify correct B-Shaped closure.

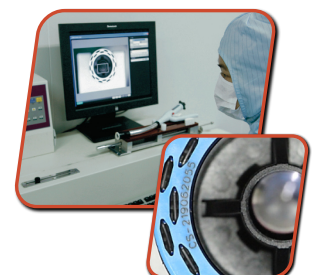

[www.frankenman.com](http://www.frankenman.com) | [info@frankenman.co.uk](mailto:info@frankenman.co.uk) | +86-512-68786150

Frankenman International Ltd - 108 Jinfeng South Road, High New District, - Suzhou, 215163, P.R. China